




CBCT Protocol

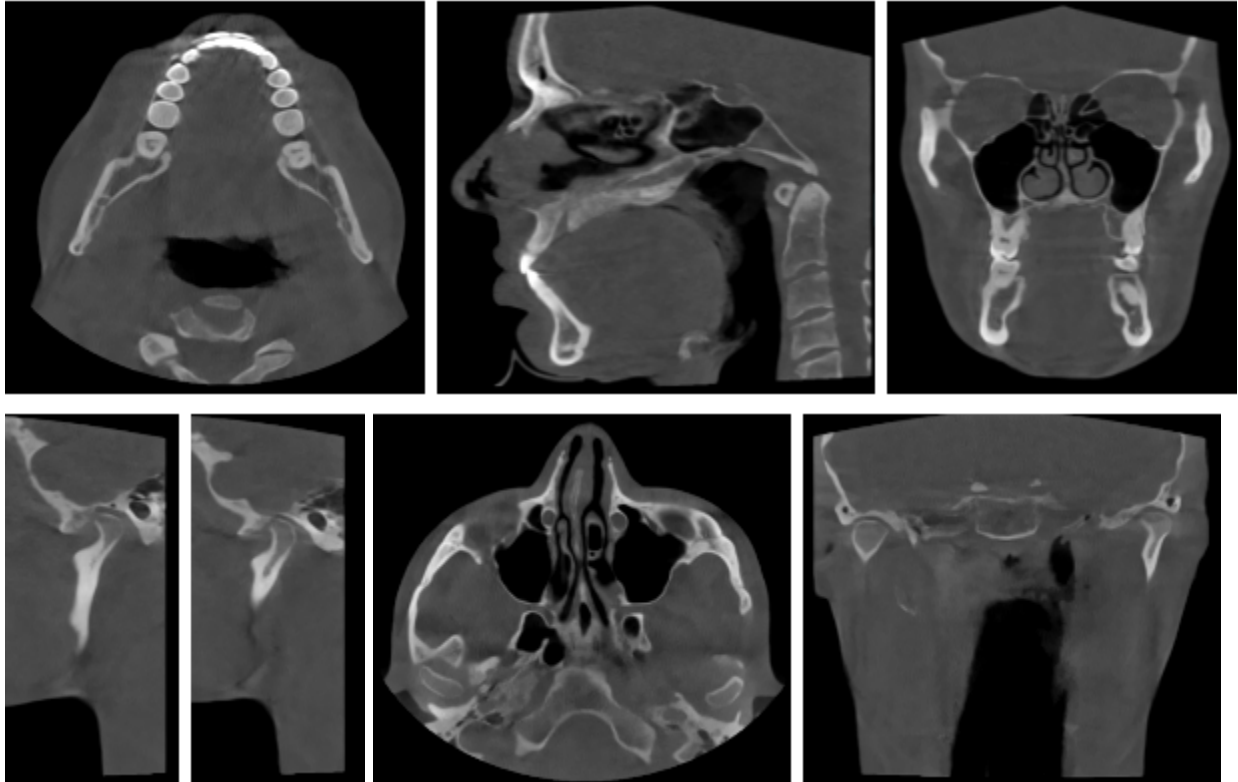

Area of Capture

- Nasion
- Tip of the nose
- TMJ
- Cervical spine (C3-C4)
- Hyoid bone

Field of View

- **Ideal:** 15x15 cm
- **Good:** 13x13 cm, 13x15 cm, 13x16 cm
- **Acceptable:** 12x12 cm, 12x13 cm, 12x15 cm, 12x16 cm

Scan Resolution

- **Voxel Size:** 0.25-0.3 mm
- **Scan Time:** Medium, 15-30 seconds

Patient Position

1. **Preparation:**
 - Remove all metal objects, including earrings, glasses, necklaces, hearing aids.
2. **Alignment:**
 - Ensure the face is in a straight, natural position, **don't bring a head forward!**
 - Instruct the patient to swallow before the scan but not during.
 - Position the teeth in habitual bite (centric occlusion or centric relation).
 - Relax the tongue.

During the Scan

- Instruct the patient to hold their teeth together.
- Encourage normal breathing throughout the scan.

This protocol ensures accurate and consistent CBCT imaging for effective diagnosis and treatment planning in sleep apnea cases.